







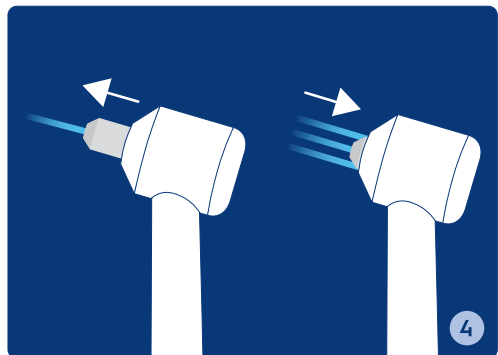
P&G RESTRICTED




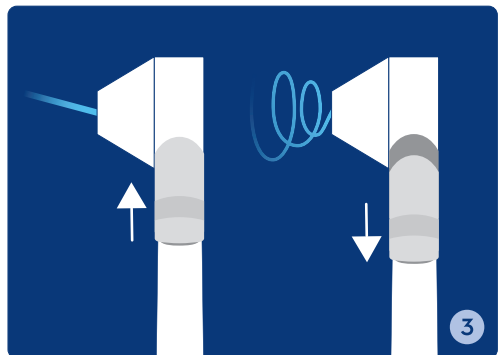
OxyJet

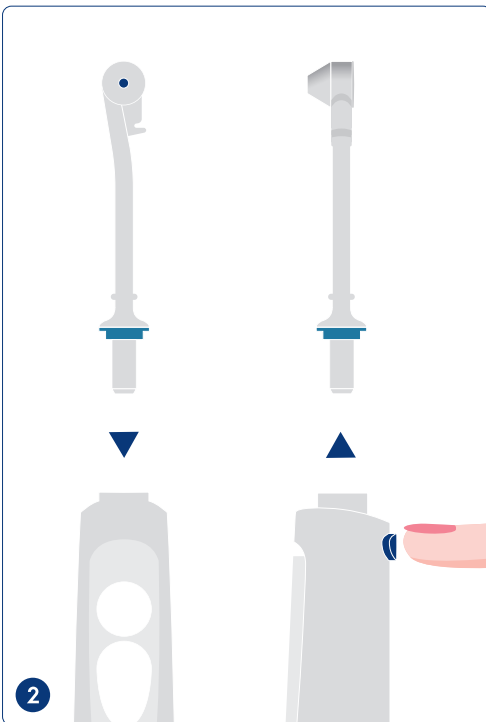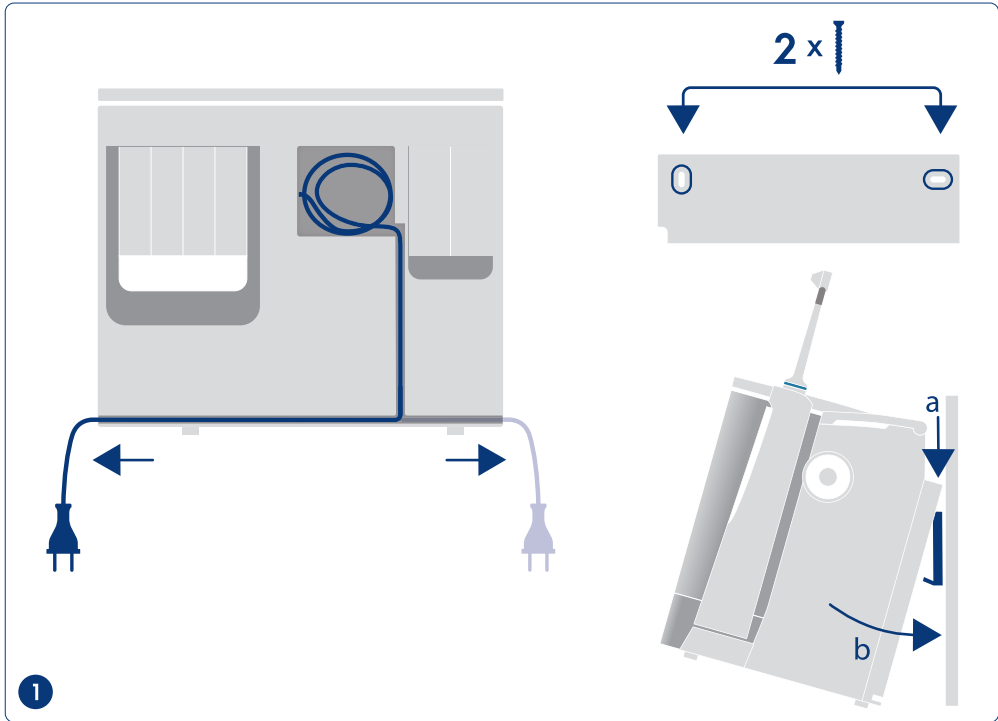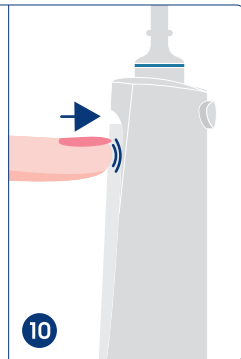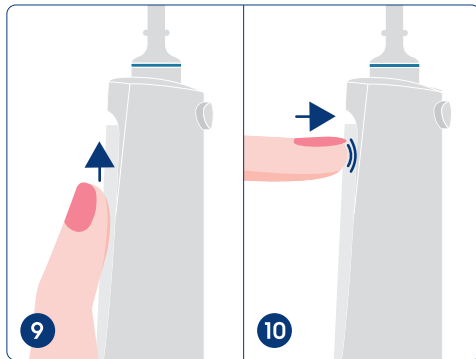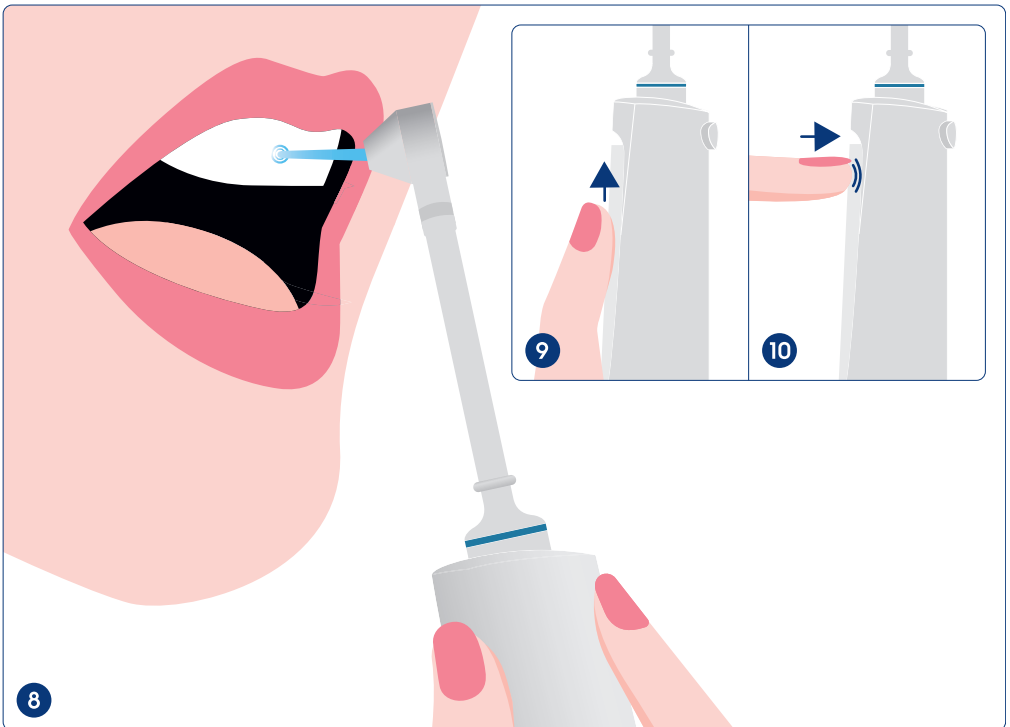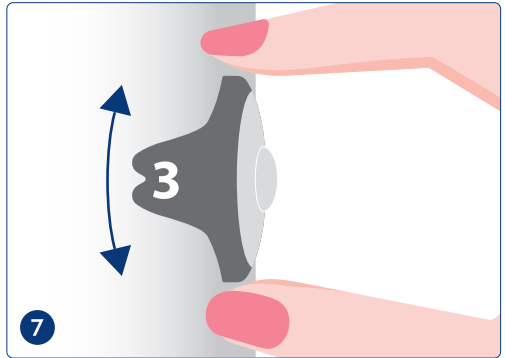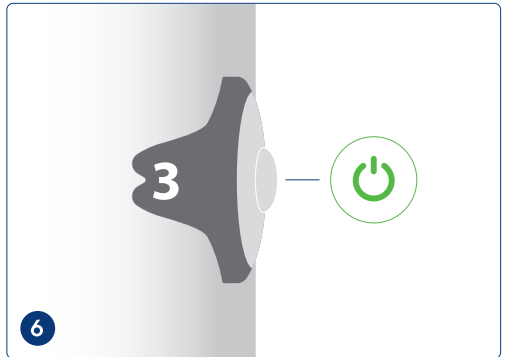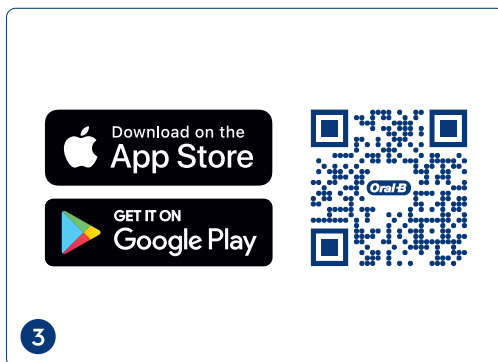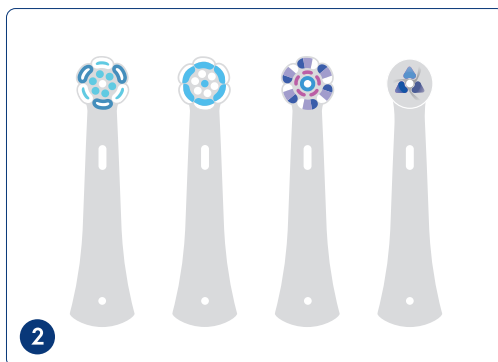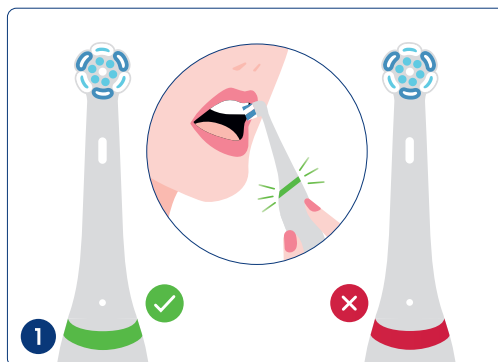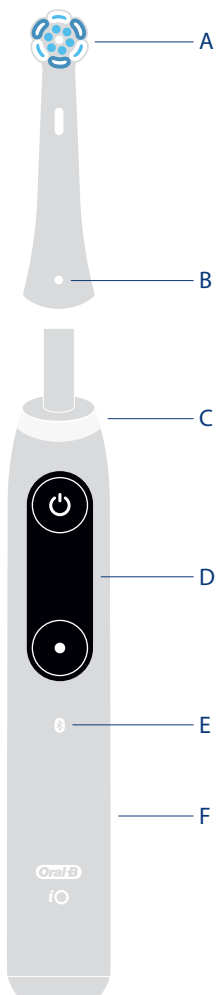Deutsch

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English

Welcome to Oral-B!

Before operating this appliance, please read these instructions and save this manual for future reference.

IMPORTANT

- The installation instructions shall state that parts that have to be fixed must be fixed so that they cannot fall into water, unless they are of IPX7 construction.
- Periodically check the entire product/cord for damage. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- If the product/ cord is damaged, take it to an Oral-B Service Centre.
- Never insert any object into any opening of the appliance.
- Oral irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- Oral irrigators shall be used only with liquids defined in the manufacturer's instructions.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in the appliance. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- The appliance must be mounted in a way that it cannot fall into water or drop to the floor. The plug should always be accessible.
- Only use charger provided with your appliance.

- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

WARNING

- Do not place the charger, smart plug or charging travel case in water or liquid or store where it can fall or be pulled into a bath or sink. Do not reach for it when fallen into water. Unplug immediately.
- Avoid sharp bending, jamming or squeezing of the cord.
- This appliance contains batteries that are non-replaceable. Do not open and disassemble device. For battery recycling please dispose whole device according to local environmental regulations. Opening the handle will destroy the appliance and invalidate the warranty.
- For tooth brushing only.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- Do not chew or bite on brush. Do not use as a toy or without adult supervision.
- Usage by children under age 3 is not recommended.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- For hygienic reasons we recommend the use of individual brush handles and brush heads for each family member.
- Small parts may come off, keep out of reach of children.
- This irrigator is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Do not chew or bite on nozzle. Do not use as a toy or without adult supervision.
- For water flossing only.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.

To avoid brush head breakage which may generate a small part choking hazard or damage teeth:

- Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fits properly. Never use without a brush head.
- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Clean brush head properly after each use (see «Cleaning Guidelines» paragraph). Proper cleaning ensures the safe use and functional life of the toothbrush.
- Replace brush (head) every 3 months or sooner if brush head becomes worn.

Important Information

- This toothbrush and pouch (accessory) contain NdFeB magnets and generate magnetic fields. If you use a medical device or have an implanted medical device, consult with the medical device manufacturer or a physician, if you have questions about safely using your Oral-B iO toothbrush.
- Your toothbrush can be used with your smartphone (find details under «Connecting your toothbrush to your Smartphone»). To avoid electromagnetic interference and/or compatibility conflicts, disable *Bluetooth* on your toothbrush handle (6) before using it in restricted surroundings, such as aircrafts or specially marked areas in hospitals.
- Persons with pacemakers should always keep the toothbrush more than 15 centimeters (six inches) from the pacemaker while turned on. Anytime you suspect interference is taking place disable *Bluetooth*.

Description

Irrigator

- A Brush head compartment
- B Toothbrush charger
- C Water container (600 ml)
- D Nozzle switch
- E Nozzle
- F Handle switch (water stop and go)
- G Irrigator handle
- H Nozzle compartment
- I On/off switch
- J Water intensity control dial

Toothbrush

- A Brush head
- B Brush head identification symbol
- C Smart Pressure Sensor / Light ring
- D Interactive display
- E *Bluetooth*®
- F Toothbrush handle
- G On/off button
- H Mode button

Specifications

For voltage specifications refer to the bottom of the charging unit.

Noise level toothbrush: 64 dB (A)

Noise level Oral irrigator: 76 dB (A)

A - Using your Oral irrigator

Your Oral irrigator comes with a convenient wall holder. It can be mounted easily using the attached items (picture 1).

Before using for the first time, fill the water container (C) and allow the water to run through completely by holding the irrigator handle (G) over a washbasin.

1. Attach one of the irrigator nozzles (E) to the handle until it locks into place (picture 2). Fill the water container with luke warm water (picture 5).
2. Turn on the irrigator by pressing the on/off switch (I) (picture 6). (When doing so, the handle switch should be in the lower position = water stop).
3. Select the water pressure with the pressure control dial (J) (1 = gentle, 5 = strong) (picture 7). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the Oral irrigator.
4. Bend over a washbasin and guide the nozzle to your teeth and gums (picture 8). Slide the handle switch (F) up to start the water stream (picture 9). Keep your mouth slightly open to allow the water to run into the washbasin.

Slide the handle switch up (F) for a continuous water stream (picture 9).

Slide the handle switch down and press into it, to use the on-demand function (picture 10).

Caution: Turn off the Oral irrigator by switching off the on/off switch once the water container is empty. Longer operation of the Oral irrigator without water may damage the appliance.

OxyJet nozzle


Nozzle functions:

You can choose between two different nozzle functions by sliding the nozzle switch (D) down or up (picture 3). Before doing so, turn off the water stream by sliding the handle switch down:

Rotating stream for general cleaning and gum massage (nozzle switch in lower position): The rotation of the water stream and its micro air-bubbles are designed for easy penetration into the interdental space and gentle gum massage.

Focused stream type (nozzle switch in upper position): The focused water stream with its micro-bubbles is designed for targeted cleaning.

WaterJet nozzle

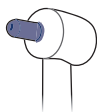Jet functions:

You can choose between different stream types by clicking-in and clicking-out the blue centre of the nozzle (picture 4). Before doing so, turn off the water stream by sliding the handle switch down:

When the blue centre is out (default position) you will get a high pressure single precision stream. As soon as you click the blue centre in, you will get a multi-jet stream, designed for great reach and massage effect of the gums and the gumline.

Maximal operation time: 15 minutes
Cooling off time: 2 hours

After use

Always empty the water container entirely. Keep the Oral irrigator switched on until the water stops running through the nozzle in order to prevent accumulation of bacteria. As soon as no more water flows out of the nozzle, switch the device off (picture 6). Press the nozzle release button (picture 2) to take off the nozzle from the Oral irrigator. Always wipe dry the handle and the Oral Irrigator; place the handle back onto its base and store the Oral irrigator in the nozzle compartment (H).

Cleaning recommendation

After use, rinse the nozzle under running water. Regularly clean the housing with a damp cloth. The water container and the nozzle compartment are removable and dishwasher-safe.

B - Using your toothbrush

First use

Please wake up your iO toothbrush by pressing the on/off button (G) and follow the instructions on the interactive display (D) to get started for first-time setup.

If not further instructed on the display, you may toggle with the mode button (H) and select with the on/off button (G).

You may change the settings via the «Settings» menu (12) at any time.

Charging and Operating

This toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- To charge the battery place the toothbrush on the plugged-in charger (B).
- The charge level is indicated on the interactive display (D). Once it is fully charged the display turns off. A full charge typically takes 16 hours.
- If the battery is running low the motor reduces its speed and the Smart Pressure Sensor (C) shortly

blinks in red. Once the battery is empty, the motor stops; it will need at least a 30-minute charge for one usage.

Note: In case the battery is empty (display off), charge for at least 30 minutes.

- You can always store the toothbrush on the plugged-in charger to maintain it at full charge; a battery overcharge is prevented by the device.

Note: Store the toothbrush at room temperature for optimal battery maintenance.

Caution: This toothbrush contains a Li-Ion battery. Do not expose the toothbrush to temperatures higher than 50 °C / 122 °F.

Using your iO toothbrush

Start Brushing

Place the brush head on the toothbrush. A small gap should remain between handle and brush head. Activate the display by picking-up the toothbrush, lifting it off the plugged-in charger or pressing the on/off button (G) (display turns on). Wet brush head and apply toothpaste. To avoid splashing, guide the brush head to your teeth before switching the toothbrush on by pressing the on/off button.

Note: Display turns off after some time. You may press the on/off button to activate display again. Guide the brush head from tooth to tooth. Brush the outsides, the insides and the chewing surfaces. Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you. After brushing, the display gives you feedback based on your brushing time.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and/or gums, Oral-B recommends using the «Sensitive» or «Super Sensitive» modes (optional in combination with an Oral-B «Gentle Care» brush head (picture 2). Availability of modes is dependent on model purchased.

Toothbrush Timer (Stutter Timer)

The iO toothbrush is equipped with timer features, which help you to achieve the professionally recommended 2-minute brushing time. The timer signals with a short stuttering sound at 30-second intervals to move to the next quadrant of your mouth. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time. The timer memorizes the reached brushing time, even when the toothbrush is briefly switched off during brushing. The timer resets if pausing longer than 30 seconds or placing the toothbrush on the plugged-in charger.

Note: During the use of the Oral-B App, the stutter

timer is deactivated in the toothbrush. You may also customize the stutter timer settings via Oral-B App.

Toothbrush Light Ring (M)

Your iO toothbrush is equipped with a light ring (C) that has various functions indicated via colours:

Smart Pressure Sensor (Picture 1)

White (default) = Operating indicator (choose individual colour via Oral-B App)/ Low pressure indicator

Green = recommended pressure indicator

Red = high pressure indicator

If the Smart Pressure Sensor illuminates in white (or your pre-selected colour), we recommend applying more pressure. The green light indicates recommended brushing pressure. If too much pressure is applied, it lights up in red reminding you to reduce pressure.

During high pressure the brush head movement will be reduced (in «Daily Clean», «Intense» and «Sensitive» modes). Periodically check the operation of the Smart Pressure Sensor by pressing moderately on the brush head during use.

Note: During the use of the Oral-B App, the Smart Pressure Sensor indication may be inactive for some features.

Note: Smart Pressure Sensor is deactivated during «Tongue Clean» mode (depending on model).

Other functions

Red = Low battery indicator

Blue = *Bluetooth* pairing indicator

Orange = Error indicator (follow display messages)

Brushing modes (depending on model)

Your iO toothbrush offers different brushing modes indicated on the display:

«Daily Clean» – Standard mode for daily cleaning

«Sensitive» – Gentle, yet thorough cleaning for sensitive areas

«Gum Care» – Gentle massage of gums

«Whiten» – Polishing for occasional or everyday use

«Intense» – For an extraordinary clean feeling (operated with a higher speed)

«Super Sensitive» – Super gentle cleaning for extra sensitive areas

«Tongue Clean» – Tongue cleaning for occasional or everyday use

Note: The sequence and availability of brushing modes may vary depending on model purchased.

When operating the «Tongue Clean» mode, we recommend using a «Gentle Care» brush head. You can brush your tongue either with or without toothpaste. Brush the whole area of your tongue systematically in gentle movements. The recommended brushing

time is 30 seconds; the motor stops after 30 seconds.

Mode setting:

Wake up your iO toothbrush to identify selected mode for your brushing session. To change to other modes, successively press the mode button (8) during use. The toothbrush automatically starts in the mode selected last.

Note: You may also customise mode settings via Oral-B App.

Interactive Display

Activate the interactive display (D) by picking up the toothbrush, lifting it off the charger or pressing the on/off button (G).

(5-11) «Mode selection» – toggle by pressing the mode button during use.

(12) «Settings» – toggle through all modes until you reach the option «Settings»; to enter press the on/off button:

- (14) manage *Bluetooth* settings (enable, disable)
- (15) select/change colour of light ring
- (16) select/change language
- (17) perform factory reset

(19) «Timer» – displays brushing time.

(20) «Battery status» – shows battery status after use, when charging or when battery is low.

(21) «Exchange reminder» – reminds you when to replace your brush head.

(22) «Brushing feedback» – brushing feedback is indicated with a smiley face.

When the toothbrush is picked up, lifted off the plugged-in charger or when you press the on/off button the toothbrush handle indicates its readiness.

«Error messages» – «outside operating limits»: you are operating or charging outside standard temperature range.

Brush heads

Oral-B offers you a variety of different Oral-B iO brush heads that only fit your Oral-B iO toothbrush handle and can be used for tooth-by-tooth precision cleaning.

Most iO Oral-B brush heads feature INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the bristle colour will fade halfway approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the colour recedes, you may exert too much

pressure on teeth and gums. In general, we recommend replacing the brush head every 3 months.

Brush Head Replacement Reminder

The toothbrush is equipped with a «Brush Head Replacement Reminder» on the interactive display, indicating the need to exchange the brush head after approximately 3 months when used twice per day.

To reset the «Brush Head Replacement Reminder» follow the instructions on the display.

Note: You may also customize the brush replacement settings via the Oral-B App.

Connecting your iO toothbrush to your Smartphone/Oral-B App

The Oral-B App offers various features:

- Personalise your toothbrush settings (e.g. choose your preferred light ring colour, adjust stutter timer, brush mode sequence or brush replacement reminder settings)
- Track your brushing behavior (e.g. brushing time and applied pressure)
- Real-time brushing recognition to help you brush every area for the right amount of time.

The Oral-B App is available for smartphones running iOS or Android™ operating systems. Check app.oralb.com for compatible smartphones and country availability of the Oral-B App. It can be downloaded for free from App Store^(SM) / Google Play store™ or you may scan the QR code on the back page of this manual (picture 3).

- Launch the Oral-B App. It will guide you through the entire *Bluetooth* pairing procedure. You may also pair new devices at any time in the Oral-B App's settings.

Note: Oral-B App functionality is limited if *Bluetooth* is deactivated on your smartphone. (For directions please refer to your smartphone's user manual).

- You may change *Bluetooth* settings on your toothbrush via the display -> Settings -> *Bluetooth* enable/disable.
- Keep your smartphone in close proximity (within 5m distance) when using it with the toothbrush. Make sure the smartphone is safely positioned in a dry place.

Note: Your smartphone must support *Bluetooth* 4.2 (or higher) / *Bluetooth* Smart in order to pair your toothbrush.

Caution: Follow your smartphone's user manual to ensure your phone/ charger is designed for use in a bathroom environment.

Cleaning Guidelines


After each brushing, rinse your brush head under running water with the handle switched on. Remove

18

the brush head and clean the handle and brush head separately. Make sure water runs inside the top of the handle while the handle is switched on, as well as inside the brush head. Wipe all parts dry and let them air-dry before reassembling your toothbrush again.

Unplug the charger before cleaning. Clean the charger regularly with a damp cloth.

Please check individual cleaning instructions for each accessory in its «Accessories» paragraph.

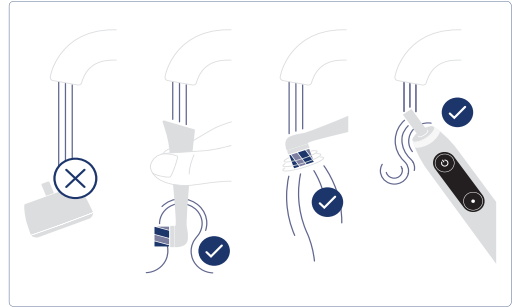Unplug the charger before cleaning. Clean the charger regularly with a damp cloth only.

Subject to change without notice.

Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase.

Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. Depending on availability, replacement of a unit might result in a different colour or an equivalent model.

This guarantee extends to every country where this appliance is supplied by Braun / Oral-B or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads/nozzles, as well as defects that have a negligible effect on the value or operation of the appliance.

The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun / Oral-B parts are not used.

These appliances are equipped with a *Bluetooth* Smart approved radio interfaces. Failure to establish a *Bluetooth* Smart connection to specific smartphones is not covered under the device guarantee, unless the malfunction of the *Bluetooth* interface

is referable to Oral-B. The *Bluetooth* devices are warranted by their own manufacturers and not by Oral-B. Oral-B does not influence or make any recommendations to device manufacturers, and therefore Oral-B does not assume any responsibility for the number of compatible devices with our *Bluetooth* interface.

Oral-B reserves the right, without any prior notice, to make any technical modifications or changes to system implementation of device features, interface, and menu structure changes that are deemed necessary to ensure that Oral-B systems function reliably. To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Customer Service Centre.

Warranty – Exclusions and Limitations

This limited warranty applies only to the new products manufactured by or for The Procter & Gamble Company, its Affiliates or Subsidiaries («P&G») that can be identified by the Braun / Oral-B trademark, trade name, or logo affixed thereto.

This limited warranty does not apply to any non-P&G product including hardware and software. P&G is not liable for any damage to or loss of any programs, data, or other information stored on any media contained within the product, or any non-P&G product or part not covered by this limited warranty.

Recovery or reinstallation of programs, data or other information is not covered under this limited warranty. This limited warranty does not apply to (i) damage caused by accident, abuse, misuse, neglect, misapplication, or non-P&G product; (ii) damage caused by service performed by anyone other than Braun or authorized Braun service location; (iii) a product or part that has been modified without the written permission of P&G.

Information

Bluetooth® radio module

Although all indicated functions on the *Bluetooth* device are supported, Oral-B does not ensure 100% reliability in connection, nor feature operation consistency.

Operation performance and connection reliability are direct consequences of each individual *Bluetooth* device, software version, as well as the operating system of those *Bluetooth* devices and company security regulations implemented on the device.

Oral-B adheres to and strictly implements the *Bluetooth* standard by which *Bluetooth* devices can communicate and function in Oral-B iO toothbrushes.

However, if device manufacturers fail to implement this standard, *Bluetooth* compatibility and features will be further compromised and the user may experience functional and feature-related issues. Please note that the software on the *Bluetooth* device can significantly affect compatibility and operation.

Warranty replacement brush heads/nozzles

Oral-B warranty will be voided if electric rechargeable toothbrush handle/irrigator damage is found to be attributed to the use of non-Oral-B replacement brush heads/nozzles.

Oral-B does not recommend the use of non-Oral-B replacement brush heads/nozzles.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads/nozzles. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads/nozzles, as communicated with the electric rechargeable handle/irrigator at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads/nozzles.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads/nozzles on handle/irrigator wear.

All Oral-B replacement brush heads/nozzles carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads/nozzles or handle/irrigator parts under any other brand name.

This guarantee in no way affects your rights under legal/statutory law.

Troubleshooting

Problem	Possible reason	Remedy
APP		
Toothbrush does not operate (properly) with Oral-B App.	1. Oral-B App is closed. 2. <i>Bluetooth</i> is not activated on smartphone. 3. <i>Bluetooth</i> is disabled on toothbrush. 4. <i>Bluetooth</i> connection to smartphone has been lost. 5. Your smartphone does not support <i>Bluetooth</i> 4.2 (or higher) / <i>Bluetooth</i> Smart. 6. Outdated Oral-B App. 7. Your smartphone is not compatible with the Oral-B App.	1. Launch Oral-B App. 2. Activate <i>Bluetooth</i> on smartphone (described in its user manual). 3. Enable <i>Bluetooth</i> via display settings on your toothbrush. 4. Keep your smartphone in close proximity when using it with the toothbrush. 5. Your smartphone must support <i>Bluetooth</i> 4.2 (or higher)/ <i>Bluetooth</i> Smart in order to pair toothbrush. 6. Download latest Oral-B App. 7. Check app.oralb.com for compatible smartphones and country availability of the Oral-B App.
Toothbrush		
Unwanted toothbrush behavior	Toothbrush was misconfigured	Perform a factory reset to restore initial settings. Navigate on toothbrush display to Settings -> Factory Reset (17). Press on/off button to confirm.
Toothbrush is not operating (during first use).	1. Battery has a very low charge level; handle does not operate. 2. Operating outside standard temperature range (message on display and/or alternating lights of light ring)	1. Charge for at least 30 minutes. 2. Always operate in standard operating temperature range (> 0°C and < 40°C).
Short stuttering sounds after 2 minutes or toothbrush does not provide any stutter timer.	1. Timer has been modified/deactivated via Oral-B App. 2. The toothbrush timer is disabled during use of brushing recognition with Oral-B App.	1. Use Oral-B App to change timer settings or perform toothbrush factory reset. 2. Disconnect toothbrush from Oral-B App.
Toothbrush does not charge.	1. Battery is at a very low charger level (0%); it may take up to 10-15min until display turns on. 2. Ambient temperature for charging is out of valid temperature range (< 0°C and > 40°C). 3. Charger may not be (properly) connected.	1. Charge for at least 30 minutes. 2. Recommended ambient temperature range for charging is 5°C to 35°C. 3. Ensure charger is connected to mains voltage.
Rattling/loud noise during brushing or turning on.	1. Brush head may not be placed correctly onto the brush handle. 2. Brush head may be worn out.	1. Check proper fit of brush head on handle. Small gap remains in between. 2. Replace brush head.

The rechargeable iO toothbrush contains a radio module which operates in the 2.4 to 2.48 GHz band with a maximum power of 1mW.

Hereby, Braun GmbH declares that the radio equipment types 3786, 3796, 3795, 3758, 3753, 3776, 3794, 3771, 3765 and 3767 are in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address: www.oralb.com/ce.

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